

TherAdapt® Products Inc.

922 N. Washington Avenue, Ludington, MI 49431

Phone: 800-261-4919 • Fax: 866-892-2478

www.theradapt.com • customerservice@theradapt.com

Platform Rocker

Please read all of the following information before using this product.

TherAdapt® Products Inc. has designed and manufactured this product with the well being of the user and caregiver(s) in mind. It is our endeavor to provide products which are attractive, comfortable, useful, highly adjustable, long lasting and safe. As when using any product, precautions must be taken.



WARNING: This product must only be used under the recommendation of a physician and / or licensed therapist and with adult supervision. Never leave user unattended. It is your responsibility to properly use and maintain your Platform Rocker. Failure to follow USE and CARE instructions or ASSEMBLE the rocker properly as outlined may result in bodily injury.

THE MANUFACTURER IS NOT RESPONSIBLE FOR INJURIES RESULTING FROM ANY CUSTOMIZATION OR MODIFICATION OF THE EQUIPMENT THAT WAS NOT COMPLETED BY THE MANUFACTURER.

USE: The Platform Rocker is designed for a single purpose: as a freestanding rocker. It is to be used by a single individual that is seated in the center of the rocker, in a manual wheelchair that is secured to the rocker as indicated in the instructions. It is to be used on the floor only. **DO NOT SUSPEND. No other use is permitted.**

CARE: Take care to inspect the ongoing condition of all parts of your Platform Rocker. Verify that all knobs, bolts, and straps are properly positioned and secure. Contact the manufacturer for replacement or repair of necessary items. All components may be cleaned with mild soap and water, or a gentle sanitizing spray. **Do not** use abrasive cleaners. Wipe wood components dry after cleaning. Excess amounts of liquid on the wood may damage the finish. Always store indoors.

WARRANTY: The manufacturer provides a full one (1) year warranty, beginning the date shipped by the manufacturer. The warranty provides full replacement of all parts against breakage when used in compliance with the manufacturer's USE and CARE instructions. Use or modification beyond the normal capacity or range of adaptability designed into the products will void the warranty. The foregoing is TherAdapt's sole liability for consequential damages.

PLEASE COMPLETE AND RETAIN FOR YOUR RECORDS

Platform Rocker PR- _____

Serial Number: _____

Purchase date: _____

***TherAdapt*[®] Platform Rocker (PR-A, PR-C, PR-E) Instructions**

Packing List:

- (1) Rocker base with (4) eye bolts
- (1) Spring loaded ramp stop and (1) slide bolt ramp stop
- (2) Wheel stops
- (4) ¼ - 20 Phillips head screws
- (2) Straps with hooks
- (2) Straps with hooks and ratcheted closure

Tools needed:

- (1) Phillips head screwdriver

The basic unit comes partially assembled. To complete assembly:

Attach the straps to the eye bolts

Place a hook from each strap through an eyebolt on one side of the Platform Rocker. Place a hook from each strap with ratcheted closure through an eyebolt on the opposite side of the Platform Rocker.



Once assembled it is ready for use-


A. Roll Wheelchair onto rocker

1. Release the rear ramp stops and lower the rear wheelchair ramp.
2. Roll a manual wheelchair or stroller frame onto the platform until the front casters strike the wheel stops. **Note:** The wheel stops can be moved forward or backward by removing the Phillips head screws and replacing them in a different set of paired “T-nuts”. **Tighten securely.**
3. Raise rear wheelchair ramp and secure with ramp stops.



B. Secure Wheelchair

1. Place one strap through the front **FRAME** of the chair and the other strap through the rear **FRAME** of the chair. **Do not** secure the chair by the footrest hanger, front casters, or rear wheels. **Do not** “loop” the straps around the chair. This will decrease your ability to tighten the straps.
2. Place the end of strap through the **bottom** of the metal buckle. Feed it forward through the cam closure (note arrow on side of buckle) and pull tight.
3. To release, push lever in center of buckle and remove strap.

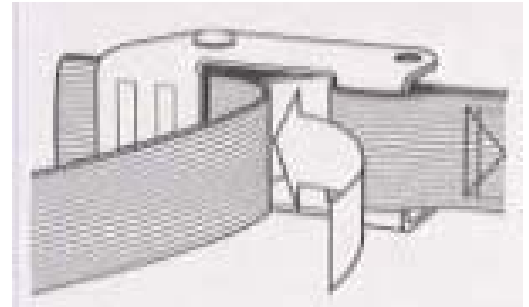

FOR PROFESSIONAL ASSISTANCE CALL 1-800-261-4919