
MODULAR CHAIRS WWW.THERADAPT.COM

# Odyssey Chair

The **Odyssey Chair** is a multifunctioning and highly adjustable chair designed to meet the needs of physically challenged children not only today but in the future as well. There are two sizes available: Early Childhood and Secondary. As an added benefit there is an upgrade package available that allows the Early Childhood Chair to grow to a Secondary size chair without purchasing a whole new chair!


The modular concept provides multiple options and adjustments with tremendous "growth" potential. It is an ideal chair for evaluation, treatment, and seating intervention in the home, school, or clinic setting. The easy adjustability of the chair can meet the changing needs of the user over the years.



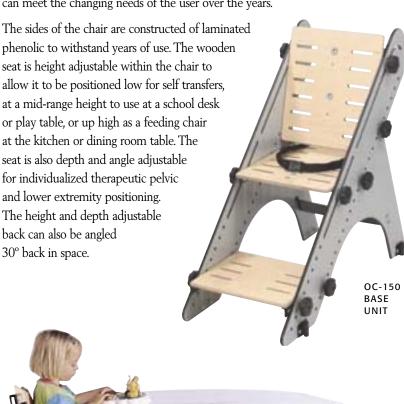

Addie, age 5, is able to play or eat at the family table in the Odyssey Chair (OC-150).



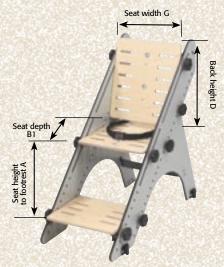

## Odyssey Chair EARLY CHILDHOOD

OC-150 Odyssey Chair

Seat depth (B1): 7" - 14"

Seat width (G): 121/2" (6" to 121/2")\*

Seat height to footrest (A): 5" to 14"
Seat height to floor: 11" - 19"\*\*

Back height (D): 15"

Overall dimensions: 29"H x 32"D x 13½W (16"W w/wheels)

User weight limit: 125 lbs.

OC-150AB Abduction Block

OC-150AR Armrests (pair)
Armrest height (C1): 6" -

OC-150BK Butterfly Kit

OC-150FS Foot Sandals (pair)

Inside dimensions: 3"W x 5"L

Outside dimensions: 41/4"W x 51/2"L

OC-150HR Headrest

Back height (E): 21" - 26"

OC-150HS Heel Stop

OC-150PS Seat and Back Pad Set

OC-150PSW Seat (full width) and Back Pad Set

OC-150SB Support Blocks\* 5"H x 4"W (pair)

OC-150SK Support Kit

(includes FS, PS, and 2 pair SB)

OC-150SW Head Support Wedges 6"H x 6"W (pair)

OC-150T Wood Tray\*\*\*

Tray dimensions: 17"W x 13"D
Tray cutout: 14"W x 4"D

OC-150TC Rounded Clear Polycarbonate Tray\*\*\*

Tray dimensions: 26"W x 21"D Tray cutout:  $13\frac{1}{2}$ "W x 6"D

OC-150WK-FL Full-Lock Wheel Kit Height: 3"H

\*Note: Seat width reduction is achieved by use of Support Blocks.

\*\*Note: Height measured without Wheel Kit.

\*\*\*Note: Armrests must be ordered to attach the Tray to the chair.

OC-UP Upgrade Kit

(Expand OC-150 to OC-250)

#### SECONDARY

OC-250 Odyssey Chair

Seat depth (B1): 11" - 16"Seat width (G):  $15^{1}/_{2}"$  (8" to  $15^{1}/_{2}"$ )\*

Seat width (G):  $15^{1}/2$ " (8" to 1 Seat height to footrest (A): 5" to 18" Seat height to floor: 11" - 19"\*\*

Back height (D): 19"

Overall dimensions: 29"H x 32"D x 17"W

(19<sup>1</sup>/<sub>2</sub>"W w/wheels)

User weight limit: 175 lbs.

OC-250AB Abduction Block

OC-250AR Armrests (pair)

Armrest height (C1): 6" - 10

OC-250BK Butterfly Kit

OC-250FS Foot Sandals (pair)

Inside dimensions: 4"W x 8"L Outside dimensions:  $5^{1/4}$ "W x  $8^{1/2}$ "L

OC-250HR Headrest

Back height (E): 25" - 30"

OC-250HS Heel Stop

OC-250PS Seat and Back Pad Set

OC-250PSW Seat (full width) and Back Pad Set

OC-250SB Support Blocks\* 5"H x 6"W (pair)

OC-250SK Support Kit

(includes FS, PS, and 2 pair SB)

OC-250SW Head Support Wedges 6"H x 6"W (pair)

OC-250T Wood Tray\*\*\*

Tray dimensions: 24"W x 22"D Tray cutout: 17"W x 4"D

OC-250TC Rounded Clear Polycarbonate Tray\*\*\*

Tray dimensions: 24"W x 22"DTray cutout:  $13^{1}/2\text{"W x }6\text{"D}$ 


OC-250WK-FL Full-Lock Wheel Kit

Height: 3"H

\*Note: Seat width reduction is achieved by use of Support Blocks.

\*\*Note: Height measured without Wheel Kit.


\*\*\*Note: Armrests must be ordered to attach the Tray to the chair.

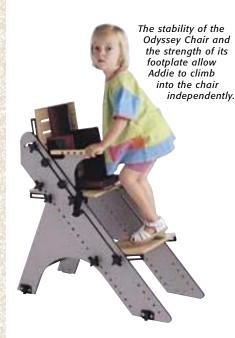







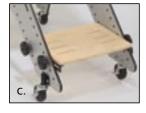

## OPTIONS

- A. The cost effective Upgrade Kit converts the Early Childhood Chair to a Secondary Chair using the components you already have. It comes complete with a new wood seat, back, and footboard, plus the hardware necessary to increase the width of the chair from 121/2" to 151/2" wide.
- **B.** A **Heel Stop**, complete dual ankle straps and curved side release buckles, is available to help keep the feet secure on the footplate. It is ideal for users who will be standing on the foot board and transferring in and out of the chair.
- C. A Full-Lock Wheel Kit, complete with mounting plates and four independent locking wheels, increases the portability of the chair from room to room, allowing the child to remain seated and well positioned for longer periods of time. The Full-Lock Wheel Kit's wheels lock both rolling and swivel movements, making it a more stable choice for independent or assisted transfers.
- D. See pages 20-21 for these options shown: Abduction Block, Armrests, Butterfly Kit, Foot Sandals, Headrest, Seat and Back Pad Set, Support Blocks (4), Head Support Wedges.



MODULAR CHAIRS WWW.THERADAPT.COM

# **Options**

# Available for: ODYSSEY CHAIRS TRANSITION CHAIRS INCLUSION CHAIRS

#### A. ABDUCTION BLOCK

A padded, adjustable depth Abduction Block is available to provide medial thigh support for children who tend to cross their legs. (Odyssey Chairs, Transition Chairs, and Inclusion Chairs)

#### **B. ARMRESTS**

Adjustable height and angle Armrests are available to provide upper extremity support and weight bearing. When ordering a tray a pair of armrests must also be ordered to attach the tray. (Odyssey Chairs, and Transition Chairs) **Note:** Armrests come standard on the Inclusion Chairs.

#### C. BUTTERFLY KIT

A Butterfly Kit is available to provide anterior chest support. It attaches with straps around the back pad of the pad set and can be adjusted in any position along the back pad for a customized fit.

Note: To be used with Pad Set or Pad Set – Wide.

#### D. FOOT SANDALS

Foot Sandals are available to secure the feet and can be placed in virtually any position needed on the footplate. They are made of molded plastic and come complete with two straps each with curved side release buckles for ease of use, comfort, and safety. (Odyssey Chairs, Transition Chairs, and Inclusion Chairs)

## E. HEADREST

A 6" high adjustable Headrest is available for additional head and neck support. It attaches to the back of the chair and has 5" of vertical adjustability. The headrest comes complete with a posterior head pad and is slotted to accept a pair of head support wedges (see J). (Odyssey Chairs, Transition Chairs, and Inclusion Chairs)

#### F. PAD SET

The Pad Set consists of a seat and back pad made of washable vinyl. The seat pad is attached with Velcro\*-like hook and loop while the back pad is secured with two bolts through the chair back. (Odyssey Chairs, Transition Chairs, and Inclusion Chairs)














G.



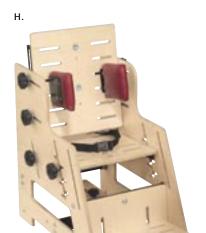

١.











#### G. PAD SET-WIDE

The Pad Set-Wide features a full width seat pad and standard back pad made of washable vinyl. The seat pad is attached with Velcro®-like hook and loop while the back pad is secured with two bolts through the chair back. Please note: Support blocks cannot be used in the seat slots when the Pad Set - Wide is in place. Support block placement in the back slots is not affected. (Odyssey Chairs, Transition Chairs, and Inclusion Chairs)

#### H. SUPPORT BLOCKS

A pair of padded, vinyl covered Support Blocks are available to provide lateral pelvic, trunk, or thigh support, or to narrow the width of the chair. They attach to the seat or back of the chair with hardware provided and can be placed in any position necessary. Several Support Blocks may be used in customizing the fit of the chair for more involved users. (Odyssey Chairs, Transition Chairs, and Inclusion Chairs)

#### I. SUPPORT KIT

A Support Kit is available at a reduced cost. The kit includes a pair of Foot Sandals, a Seat and Back Pad Set, and two pairs of Support Blocks. (Odyssey Chairs, Transition Chairs, and Inclusion Chairs)

#### J. HEAD SUPPORT WEDGES

45° angle Head Support Wedges are available to provide lateral head control. They attach to the headrest to create a cradle and can be adjusted for a customized fit. The set comes complete with a pair of vinyl pads that attach to the wooden wedge shaped mounting brackets. (Odyssey Chairs, Transition Chairs, and Inclusion Chairs)

#### K. WOOD TRAY

A rugged Wood Tray is available for functional activities and provides additional upper extremity support and positioning. The spring loaded, quick release hardware on the Odyssey Chair and Transition Chair Trays allow for easy attachment to and removal from the armrests. **Please note:** Armrests must be ordered to attach the tray to these chairs. The Inclusion Chairs come with padded armrests that allow the Tray, with "Z' shaped mounting hardware, to slide on and off.

#### L. POLYCARBONATE TRAY

A clear Polycarbonate Tray is also available. The clear tray allows the user to see his or her lower extremities for greater body awareness and body scheme. The spring loaded, quick release hardware on the Odyssey Chair and Transition Chair Trays allows for easy attachment to and removal from the armrests. **Please note:** Armrests must be ordered to attach the tray to these chairs. Also note: The polycarbonate Odyssey Chair Tray is kidney shaped and cannot be used with other chairs.