

Bolster Chair

The TherAdapt® **Bolster Chair**, which comes in three sizes, provides stability and symmetry for children whose legs tend to “cross” or adduct. The width of the flat padded “bolster bench” allows the child’s pelvis to be centered on the seating surface and blocks him from “rolling” off to one side, which can sometimes be a problem on a traditional round bolster.

The TherAdapt® Bolster Chair is height adjustable to allow for appropriate lower extremity positioning. The foot rails keep the feet and lower legs in alignment with or without the use of the standard foot straps. It comes equipped with a pelvic positioning strap, which assists in maintaining a neutral pelvic tilt. The standard height, depth, and angle adjustable tray features spring loaded locks for quick release.


BASE UNIT



The stable, symmetrical positioning provided by the (BC-200) Bolster Chair's bolster bench helps Laura, age 4, to sit and play. Note how the adjustable footplate (BC-200FP), complete with foot sandals, secures her feet in alignment.

FIT RIGHT GUIDE

See center insert page L

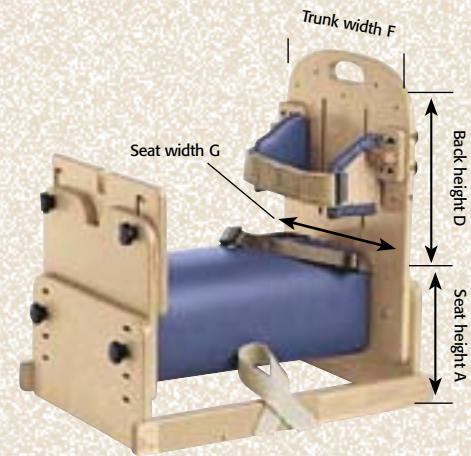

Bolster Chair

EARLY INTERVENTION

BC-100 Bolster Chair with Tray

Overall dimensions: 13"W x 22"L
 Seat width (G): 5 1/2"
 Seat height (A): 5 1/2" - 8 1/2" to floor
 Back height (D): 11"

Tray dimensions: 25 1/2"W x 24"D
 Tray cutout (F): 12"W x 4" - 7"D

BC-100TP Thoracic Pads (pair)

Trunk width (F): 4" - 6 1/2"

PRESCHOOL/PRIMARY

BC-200 Bolster Chair with Tray

Overall dimensions: 17"W x 23"L
 Seat width (G): 7"
 Seat height (A): 6 1/2" - 11" to floor
 (6 1/2" - 9 1/2" to footplate)
 Back height (D): 13"

Tray dimensions: 27 1/2"W x 24"D
 Tray cutout (F): 15"W x 5 1/2" - 8 1/2"D

BC-200FP Footplate with Sandals (pair)

BC-200TP Thoracic Pads (pair)

Trunk width (F): 6 1/2" - 10 1/2"

PRIMARY/INTERMEDIATE

BC-300 Bolster Chair with Tray


Overall dimensions: 19"W x 27"L
 Seat width (G): 8 1/2"
 Seat height (A): 9" - 15" to floor
 (9" - 13 1/2" to footplate)
 Back height (D): 17"

Tray dimensions: 29"W x 28"D
 Tray cutout (F): 17 1/2"W x 5 1/2" - 8 1/2"D

BC-300FP Footplate with Sandals (pair)

BC-300TP Thoracic Pads (pair)

Trunk width (F): 10 1/2" - 14 1/2"

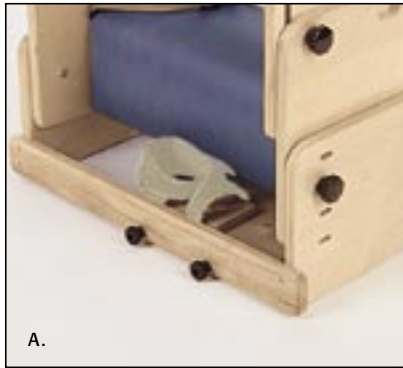

A.



B.



C.

For BC-100 order MB-A
For BC-200 order MB-A
For BC-300 order MB-B




D1.



D2.



D3.

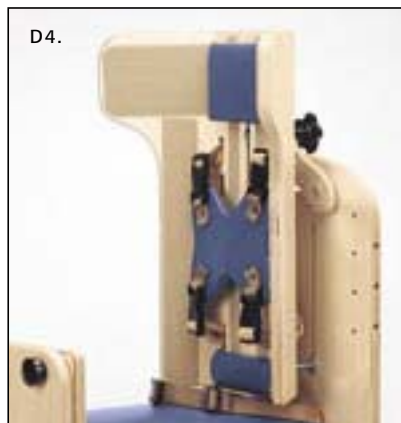

D4.

OPTIONS

- A.** For additional control a **Footplate**, complete with foot sandals, is available for the BC-200 and BC-300, to secure the lower extremities in virtually any position necessary along the base.
- B.** The **Thoracic Pads**, complete with anterior chest strap, are tapered for an intimate, comfortable fit. They can be adjusted in height and width for a variety of users.

ACCESSORIES

- C.** A **Mobile Base**, with locking caster wheels, is available to transport the Bolster Chairs. The **BC-100** interfaces with **MB-A**. The **BC-200**, with or without a Footplate, also interfaces with the **MB-A**. The **BC-300** with or without a Footplate, interfaces with the **MB-B**. Please note that the addition of a Mobile Base will raise the chair 3 1/2". For more information see page 93.
- D.** The paired slots in the Bolster Chair's back accept either the **Corner Back Insert**, the **Low Back Insert**, the **Flat Back Insert**, or the **Winged Back Insert**. These customizing options offer head and/or trunk support for the more involved user.
- D1.** Corner Back Insert. For information see page 12.
- D2.** Low Back Insert. For information see page 12.
- D3.** Flat Back Insert. For information see page 10.
- D4.** Winged Back Insert. For information see page 11.