
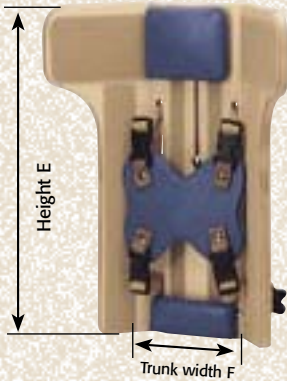


# Winged Back Insert

## FIT RIGHT GUIDE

See center insert page L

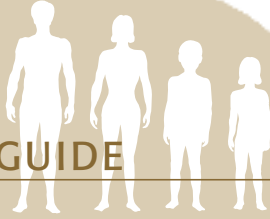

BASE UNIT

### Winged Back Insert

WBI-100 – Early Intervention

Height (E): 16" - 24"  
Trunk Width (F): 9"

WBI-100FH Flat Head Pads (pair)

WBI-100PP Protraction/Pelvic Pads (pair)

WBI-100TH Triangle Head Pads (pair)

WBI-100TT Triangle Trunk Pads (pair)

BLP-102 2" high Back Liner Pad

BLP-104 4" high Back Liner Pad

WBI-200 – Preschool/Primary

Height (E): 19" - 29"  
Trunk Width (F): 11"

WBI-200FH Flat Head Pads (pair)

WBI-200PP Protraction/Pelvic Pads (pair)

WBI-200TH Triangle Head Pads (pair)

WBI-200TT Triangle Trunk Pads (pair)

BLP-202 2" high Back Liner Pad

BLP-204 4" high Back Liner Pad

WBI-300 – Primary/Intermediate

Height (E): 23" - 35"  
Trunk Width (F): 14"

WBI-300FH Flat Head Pads (pair)

WBI-300PP Protraction/Pelvic Pads (pair)

WBI-300TH Triangle Head Pads (pair)

WBI-300TT Triangle Trunk Pads (pair)

BLP-303 3" high Back Liner Pad

BLP-305 5" high Back Liner Pad

The TherAdapt® **Winged Back Insert**, available in three sizes, was designed to interface with the Adjustable Positioning Chairs and Bolster Chairs, allowing for seat depth reduction as needed. This Insert is effective at providing thoracic and head support for the mild, moderate, or severely involved user who typically hyperextends at the trunk and retracts the scapulae. With the pelvis properly positioned in the appropriate chair, the Winged Back Insert encourages scapular protraction and normal spinal alignment.

The Insert is both height and depth adjustable. It comes standard with a mounting yoke, butterfly anterior chest support, a posterior head pad, and lumbar pad.



A.



B.



C.



D.



E.



### OPTIONS

**A. Back Liner Pads** are available in two heights to line the wood down the center of the Insert for comfort and pressure relief over bony prominences. To completely line the back of the Insert, order two additional smaller and two additional larger **Back Liner Pads**.

**B. Flat Head Pads** are available to provide a cushioned "cradle" for the winged portion of the head support.

**C. Triangle Head Pads** are recommended when lateral head support is needed. They are ideal for securing the head to promote visual gaze and tracking.

**D. The Triangle Trunk Pads** are adjustable both vertically and horizontally to provide lateral trunk support and stabilization at 90° to the trunk.

**E. Protraction/Pelvic Pads** are 25° angle pads that can be used behind the shoulder blades to encourage protraction of the scapulae, or at the pelvis for lateral positioning and stabilization, especially with a narrow pelvis.